

Xenon Sampling System for LUX

A. Dobi^a

^a*Physics Department, University of Maryland, College Park MD, USA*

Abstract

We describe a xenon purity analysis system we have developed and used for the LUX dark matter experiment based on a mass spectrometry technique. The device is fully automated, simple, compact and is integrated into the LUX circulation system allowing for hourly, in-situ sampling from several ports. The sensitivity of the spectrometer is enhanced by several orders of magnitude by the presence of a liquid nitrogen cold trap, and many impurity species of interest can be detected at the level of one part-per-billion or better. In the case of Kr, a troublesome internal background for xenon based dark matter experiments, the sensitivity is sub one part-per-trillion. We have used the technique to screen the LUX xenon before, during, and after the first underground science run, and these measurements have proven useful. This is the second application of the cold trap mass spectrometry technique to an operating physics experiment, first employed for EXO-200.

Keywords:

Noble gas, xenon, purification, SAES purifier, cold trap, mass spectrometry, mass spectroscopy, Dark Matter

Contents

1	Introduction	3
2	Purity Requirements	4
3	Methodology	4
4	Purity Figure of Merit and Calibration	8
4.1	RGA Partial Pressure Measurement	9
4.2	Correcting for RGA Gain Drift	9
4.3	Correcting for Flow Rate	11
4.4	Gain and Flow Corrected Figure of Merit	12
4.5	O ₂ and H ₂ Figure of merit	14
5	Applications of Gas Purity Analysis	16
5.1	Xenon Inventory Screening	17
5.2	⁸⁵ Kr Monitoring	17

1. Introduction

The LUX collaboration has conducted the first dark matter search at the Sanford Underground research facility in South Dakota becoming the leading spin-independent WIMP cross section limit [ref]. The LUX detector is a dual phase xenon TPC in which the energy deposition and the 3D coordinate of an event can be reconstructed from the primary scintillation (S1) and the subsequent secondary scintillation (S2) collected on two arrays with 61 PMTs each [LUX detector paper]. To ensure that radioactive backgrounds in the detector materials do not obscure the WIMP signal, a comprehensive materials screening program was employed during the construction of the experiment [BG paper]. This program measured and certified the radiopurity of all passive detector materials, including the PMTs, TPC instrumentation, the cabling, the xenon vessel, and the shielding materials. In this article we describe a complementary program to study the purity of the xenon source itself. This campaign has been carried out with in-situ measurements during the surface run in 2012 and the underground science run in 2013. These measurements have allowed us to verify that the xenon stockpiles after gas chromatography [ref] are suitable for their intended purposes in LUX, to monitor the performance of the xenon gas purifiers, to measure emission rates from plastic components, and to independently monitor the radioactive backgrounds from ^{85}Kr and ^{39}Ar beta decay daily.[ref]

The measurements described in this article are based on a cold trap mass spectrometry technique which was developed to study O_2 , N_2 , and CH_4 impurities in xenon [CT, Purifier, tritium paper]. The method has been extended for the purposes of krypton detection as well [6] and was first implemented for xenon screening use in EXO-200. [EXO paper]

Mass spectrometry has several attractive features as an analysis technique. It can detect both electronegative and non-electronegative impurity species, including O_2 and some of the problematic noble gases which contain radioactive isotopes such as Kr. It allows each impurity species to be identified and counted individually and simultaneously unlike mass spec techniques that can only measure one species at a time[ref]. It also allows the purity of the xenon to be screened, and possibly corrected, prior to detector operations, which is similar to how all other EXO-200 and LUX detector materials are treated. Further, the technique is simple and cost effective only requiring readily available laboratory equipment. The analysis system takes up minimal space and can be plumbed in directly to any circulation panel allowing for in-situ, hourly purity results at the sub $\cdot 10^{-9}$ g/g level for multiple impurity species, making the method attractive over other more complex techniques which only detect a single species at a time, off site and on a significantly larger time scale [Aprile, and the Germans].

We have found many useful applications for having an in-situ automated xenon sampling device which will be described in detail. Mainly for LUX, we could ensure the Kr content of the xenon each time the xenon was cryo-pumped, compressed or any other time when there was a potential for air contamination. Xenon100 suffered a year long setback due to a minor leak in the reticulation pump which took days to show up in the liquid TPC via data analysis, something that would have been

caught by the daily gas sampling program developed LUX. Just one liter of air breaching the plumbing is enough to exceed the radioactive background goals for the LUX experiment. It is important to note that minimal xenon is lost during the sampling which consumes only milligrams of xenon per sample, the bulk of the sampled xenon is simply recovered after analysis.

2. Purity Requirements

There are two aspects of xenon purity requirements that can be screened for using the gas sampling technique. First, the xenon stockpiles should be relatively free of electronegative impurities for successful TPC operations. For example, the O_2 concentration should be less than 1 ppb in order to collect ionization over the 20 cm drift length of the detector[8]. Other electronegative impurity species may have larger or smaller attachment coefficients. For example, the attachment coefficient of N_2 , measured in liquid Ar, is about a factor of 1000 smaller than that of O_2 [9]. Second, Krypton and argon impurities are problematic because they include the beta emitters ^{85}Kr and ^{39}Ar . After gas chromatography, LUX requires a Kr and Ar concentration of $<5 \cdot 10^{-12}$ g/g and $<1 \cdot 10^{-6}$ g/g respectively. With Q values of 687 keV and 565 keV respectively, these beta decays occur in the fiducial volume defeating the power of xenon self-shielding and contribute to the WIMP search background. Clearly the krypton goal is much more stringent than the argon goal, so krypton contamination is the more serious concern. Note that the mass spectrometry technique does not detect ^{85}Kr and ^{39}Ar directly, but instead detects the stable and abundant isotopes ^{84}Kr , ^{86}Kr , ^{82}Kr and ^{40}Ar , using the standard isotopic abundances ($\sim 2 \cdot 10^{-11}$ $^{85}Kr/natKr$ [10] and $\sim 8 \cdot 10^{-16}$ $^{39}Ar/natAr$ [11]). The concentrations of the radioactive components can then be inferred from the isotopic abundances, when these are known. Radon is another radioactive noble gas which is a serious concern. However, radon cannot be detected by our method, but a limit can be place on external radon diffusion by knowing the activity in the surrounding air and measuring Kr and Ar accumulation in the xenon overtime.

3. Methodology

The xenon gas analysis technique used for LUX is similar to the method described in Refs [5] [6] [EXO]. We use a residual gas analyzer (SRS RGA-200) mass spectrometer aided by a liquid nitrogen coldtrap to measure the impurity content of xenon gas, the diagram representing the basic setup is show in figure 1. The RGA requires a high vacuum to operate (10^{-5} torr or less), so only modest amounts of xenon gas can be admitted into the device at one time. This is accomplished with a precise vacuum leak valve. The RGA measures the partial pressures of the relevant atomic masses, and these partial pressures are proportional to both the concentration of the various gas species and to their flow rate through the leak valve. By opening the leak valve further and further, we can increase the flow rate to an arbitrarily high value, resulting in higher and higher RGA partial pressures for all species. In principle this allows very small concentrations of impurities to be detected

above background levels. In practice, however, the use of very high flow rates will cause the RGA to saturate due to the high partial pressure of the bulk xenon gas. Once the pressure of the xenon rises above 10^{-5} torr, the RGA will be unable to operate, and this limits the impurity sensitivity of the RGA to about one part-per-million. Since we desire to detect impurities at the part-per-billion level or better, we must prevent RGA saturation by removing the bulk xenon from the sampled gas. This is accomplished by placing a liquid nitrogen cold trap between the leak valve and the RGA. In the cold trap, the dominant xenon partial pressure is held fixed by the liquid nitrogen bath at the vapor pressure of xenon ice at 77 K (1.8×10^3 torr), independent of the leak valve setting and the flow rate. With these conditions in the coldtrap impurity species of interest such as N_2 , O_2 , He, Ar, Kr, CH_4 remain gaseous and pass through unabated while the bulk xenon is captured [ref]. The xenon partial pressure is further reduced below 10^{-5} torr at the RGA by including vacuum plumbing with impedance after the cold trap. With the xenon pressure now held fixed, we can increase the flow rate, and therefore the partial pressures of many impurity species, by several orders of magnitude. Depending on flow, the xenon leaving the coldtrap has up to six orders of magnitude enhancement in the impurity to xenon ratio which is then passed to the RGA for measurement. In laboratory bench tests this method has achieved sensitivity to 0.12 ppb of O_2 , 0.3 ppt of Kr (g/g) and 5.0 ppt of CH_4 . Impurities which do not pass through the cold trap at liquid nitrogen temperature, such as H_2O or heavy hydrocarbons, cannot be detected in this way, but in some cases they may be detectable as an excess background level in the the cold trap plumbing after it returns to room temperature.

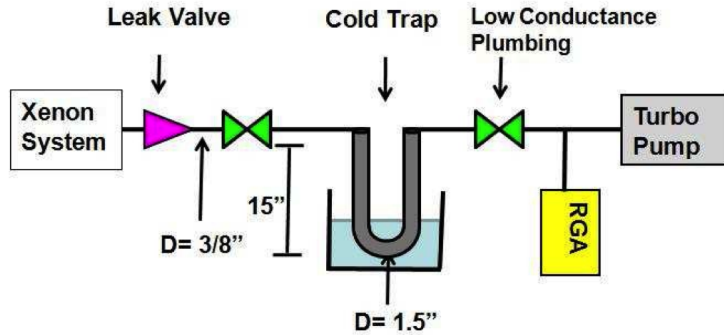


Figure 1: Basic setup of the purity analysis technique. A coldtrap is used to remove the bulk xenon, enhancing the impurity to xenon ratio of several species of interest by up to six orders of magnitude before passing the sample to a RGA for analysis.

For the LUX experiment we have integrated the xenon gas analysis system into the xenon circulation panel and automated the sampling procedure using pneumatic value via slow control, the plumbing schematic is shown in Figure 2. Xenon gas is collected from the various locations of interest along the circulation path including the input/output of the xenon gas purifier, the output of the

TPC, directly from the bottle farm and the xenon storage vessel. The gas sampling is done directly from the port of interest with the lines permanently plumbed in, reducing potential for air leaks and other systematics associated with having to draw samples into bottles and shipping them elsewhere for analysis [ref]. In the case where we sample a high pressure source such as a gas cylinder or SRV (storage recovery vessel), the sampling port is located at the output of a pressure regulator. The pressure of the xenon gas sampled is set either by the regulator (for a high pressure source) or by the pressure of the LUX xenon gas system which does not exceed 3 atm absolute. In either case the sample pressure is typically between 1 and 1.5 atmospheres in the 0.80 liter sample volume (about 5 to 9 grams of xenon). A diagram showing sampling ports plumbed into the LUX experiment and a schematic of the analysis system is shown in Figure 2.

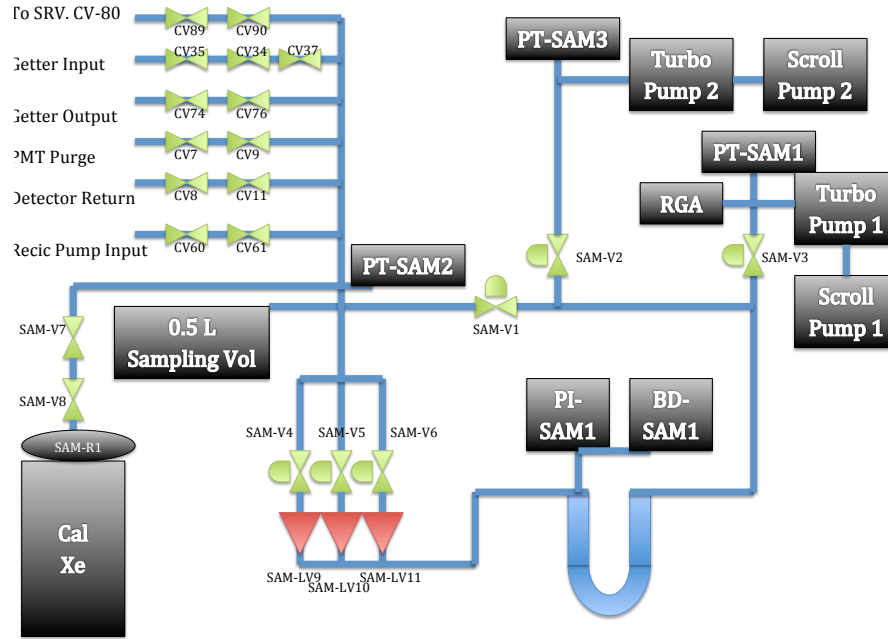


Figure 2: A schematic of the purity analysis system, the total sampling volume is 0.80 liters. Typical sample pressures are between 1 and 1.5 atmospheres.

The analysis system consists of a 0.5 l bottle used as a buffer volume for sampled xenon and 0.3 l of additional plumbing to the sampling ports, a capacitive manometer, three vacuum leak valves set at various leak rates, a U- shaped liquid nitrogen cold trap, a low-conductance plumbing element, an SRS RGA-200, two cold cathode vacuum gauges, two turbo pumps, two scroll pumps, six pneumatic VCR valves and instrumentation necessary for automation via a slow control interface. The capacitance

manometer (MKS Baratron model 627B) measures the pressure in the sample buffer volume (0.8 l total) before and during the impurity measurement. The leak valves (part number of the ones with the dial?) admit the sampled gas to the cold trap at much reduced pressure and allows for control of the flow rate into the trap at three preset values. The cold trap consists of a U-shaped segment of 3.8 cm OD stainless steel tube with a height of 38 cm designed to be partially immersed in a liquid nitrogen bath. The low conductance element is a 10 cm segment of 0.6 cm and 0.95 cm OD plumbing with one right angle bends. A port leading to the SRV is also plumbed in to the analysis system allowing for for 99.9% of the 5 to 9 grams of sampled xenon to be recovered after the conclusion of the impurity measurement, only several milligrams of xenon are consumed for a measurements. The cold cathode gauge allows the total pressure at the RGA to be monitored and recorded. We have found it necessary to limit this pressure to less than 1.0×10^{-5} torr to prevent RGA saturation effects. The analysis system includes a port to allow for a sample bottle filled with calibration xenon to be attached. We prepared the calibration xenon gas with known concentrations of O₂, N₂, He, Ar, Kr, and CH₄, by mixing known quantities of these impurities with a known amount of xenon. Except where otherwise noted, all of our impurity measurements are quantified by comparing to such calibration xenon. Once a xenon gas sample has been collected we analyze the sample for impurities as follows. The RGA is prepared to measure and record the partial pressures as a function of time of the relevant species, an RGA scan of a calibration is shown in Figure 3. We typically monitor, in atomic mass units, 28 (N₂), 32 (O₂), 18 (H₂O), 132 (Xe), 2 (H₂), 4 (He), 15 (CH₄), 40 (Ar), and 82, 84, 86 (Kr). We cool the cold trap with liquid nitrogen while it is still being pumped to ultra-high vacuum with the turbo pump, and then the xenon gas is admitted into the cold trap in three steps. First, the leak valve at the lowest flow setting is opened which admits xenon into the cold trap at a very small flow rate, below 10^{-4} standard liters per minute (SLPM). Xenon ice forms in the cold trap and the analysis system is purged of residual trace impurities by the flowing gas. We wait ten minutes for the backgrounds to stabilize before proceeding. Second, we open the leak valve to a larger flow rate, usually corresponding to an flow rate of ~ 0.01 SLPM. This coarse measurement is performed to ensure that no N₂ from an air leak or other impurity is present at amounts that could saturate the RGA signal, if the xenon is deemed pure then the system moves on to the next higher flow stage for the $1 \cdot 10^{-12}$ g/g Kr measurement. Third, the main leak valve is opened corresponding to flow rates of ~ 0.1 SLPM gaining another factor of 10x in sensitivity from the coarse measurement. The higher flow rate used in this step allows impurities to be observed with high sensitivity, typically one part-per-billion to one part-per-trillion or better depending on the impurity species. The xenon gas in the sample bottle is slowly depleted as it flows into the cold trap. This causes the sample bottle pressure, and the flow rate through the leak valve, to decrease in time while the measurement takes place (see Figure 3). After a measurement period of three minutes, the leak valve is closed, the cold trap is allowed to warm to room temperature, and the xenon in the cold trap and the sample bottle is collected in the recovery bottle with liquid nitrogen. After the sample we take an RGA mass spectrum scan (1-200 amu) of the cold trap at room temperature to look for impurity species such

as H_2O , solvents, or heavy hydrocarbons which would be trapped in the cold trap at liquid nitrogen temperature. Note that to achieve the best sensitivity, the highest possible flow rate should be used for each measurement, limited only by RGA saturation effects.

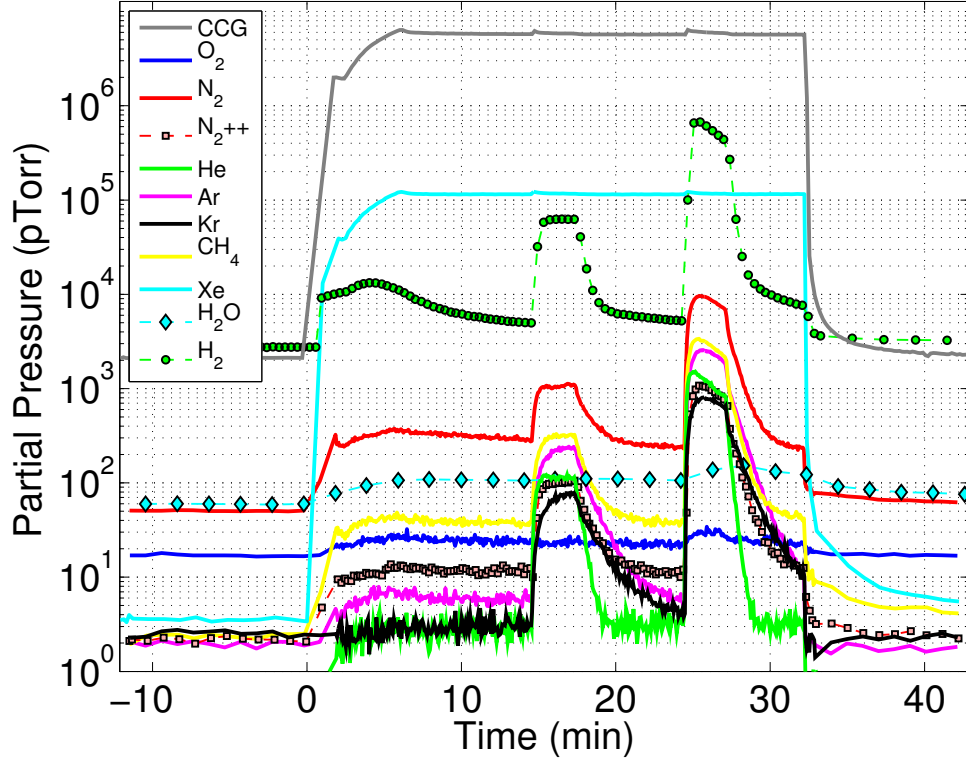


Figure 3: Typical RGA calibration measurement. (CCG pressure appears higher because of a gas correction factor of 3x for xenon, 132 amu is only 1/4 of the total naturally occurring xenon and because half of xenon is double ionized appearing at 66 amu).

4. Purity Figure of Merit and Calibration

For a given leak rate the measured RGA partial pressures have been found to be proportional to the impurity concentration over several orders of magnitude [ref]. In practice samples are taken at a variety of initial pressures and the RGA gain is prone to small, daily gain drifts. To account for the variation in initial sample pressure, which dictates flow rate through the pre set leak valve, the partial pressure signals are normalized to the average flow during the time that measurement is made. The reason being that the signal at a given atomic mass depends linearly on the number of impurities

passing by per unit time. We have shown in pervious calibrations that the leak rate normalized RGA partial pressure is proportional to the impurity concentration of the xenon gas, both in the case of constant flow or declining flow [ref ref]. We also find for the LUX purity analysis system, that the leak rate normalized partial pressure is well behaved and is discussed further in this section. Due to power outages and environmental variables the RGA is subject to daily gain shifts. The gain can be corrected for by simply calibrating the device each time the RGA is power cycled. However, frequent calibrations may not be practical over months of running so we present a convenient method to normalize the RGA’s gain by tracking xenon vapor pressure.

4.1. RGA Partial Pressure Measurement

A typical quadrupole mass spectrometer (such as the SRS RGA-200) measures the partial pressures of impurity species under vacuum conditions by trapping ionized particles at a given frequency for which the charge to mass ratio is in resonance. The particles are ionized by electrons released through thermionic emission from a hot filament and then accelerated through a 70-90 V potential. Due to the bombardment of electrons, atoms and diatomic molecules, such as Xe and N₂, may become doubly ionized and appear at half of the expected charge to mass ratio. One the other hand, molecules such as methane and carbon monoxide can be cracked and show up at masses corresponding to a lost hydrogen or split molecule, respectively. For the measurements reported in this work we track the masses given in Table 1. For a select species we report on the typical amount that a diatomic molecule or atom is doubly ionized by the RGA at 90eV in Table 2. Knowing these parameters is crucial for ensuring that measurements for the expected N₂ signal at 28 AMU is really nitrogen and not carbon monoxide or that the expected methane signal at 16 AMU is not contaminated by O₂⁺⁺. For the case of methane we track 15 AMU corresponding to CH₃⁺, which makes up a sufficient 40% of the total methane signature without the potential for contamination at 16 amu from O₂ or cracked H₂O. Xenon and Krypton, being relatively easily ionizable, show up 40% and 20% doubly ionized, respectively. For the most important specie of interst, Kr, we measure not only the most abundant isotope at 84 amu (57.0%) but also the next two significant isotopes 86 amu (17.3%) and 82 amu (11.6%). Measuring 82+84+86 instead of simply 84 was found to boost the Kr sensitivity by exactly what was expected by counting the additional isotopes, an increase of 50%. Additionally, using the ratios of the krypton isotopes is a powerful tool for discriminating potential hydrocarbon backgrounds from a legitimate krypton signal. Though not implemented in this work, another 20% gain in krypton sensitivity can be achieved by looking for the doubly ionized Kr⁸⁴ signal at 42 amu and Kr⁸⁶ 43 amu, doubly ionized Kr⁸² at 41 amu may have contamination from the argon peak at 40 amu.

4.2. Correcting for RGA Gain Drift

Before discussing calibration and purity measurements we must consider the stability of the RGA signals which can be effected by environmental variables. To our advantage there is a convenient way to monitor the gain of the RGA by tracking the one constant linking all xenon samples, the xenon

AMU	Primary Use	Secondary Use
2	H ₂	HydroCarbon, H ₂ O
4	He	
12	C	CO, HydroCarbon
14	N ₂ ⁺⁺	
15	methane (40%)	methane (45%), H ₂ O, CO
16	O ₂ ⁺⁺	
18	H ₂ O	
28	N ₂	CO
32	O ₂	
40	Ar	
44	CO ₂	
82	Kr (11.6%)	
84	Kr (57.0%)	
86	Kr (17.3 %)	
132	Xe (26.9%)	

Table 1: Partial pressures of impurity species (in AMU) tracked by the RGA during the purity analysis.

Species	AMU	++ AMU	++ Fraction
N ₂	28	14	10%
O ₂	32	16	10%
Ar	40	20	10%
Kr	84 (57.0%)	42	20%
Xe	132 (26.9%)	66	40%

Table 2: Typical values for the fraction of doubly ionized atoms and diatomic molecules at the RGA’s electron acceleration potential of 90 V.

vapor pressure. Since a liquid nitrogen bath is used for the impurity detection technique the dominant xenon pressure is always fixed in the coldtrap at 1.8 mtorr and 77K [ref]. The xenon vapor pressure is independently checked with a cold cathode gauge near the turbo pump, and is found to be stable within 3% over the course of measurements from January 2013 to January 2014. Thus, any changes in the xenon RGA measurement while pumping on xenon ice (as long as the CCG and turbo pump current are fixed) is an indication of a shift in the RGA’s gain. We corrected for daily changes in the RGA gain by normalizing the RGA signals to the xenon pressure of the measurement. Figure 4

shows the average xenon pressure measured by the RGA for each measurement and calibration over the course of a year. There are variations up to a factor of ten in the RGA's response that need to be corrected for, it is not sufficient to simply use the nearest calibration point. Figure 6 shows the difference normalizing the measurement to xenon vapor pressure makes in the signal stability.

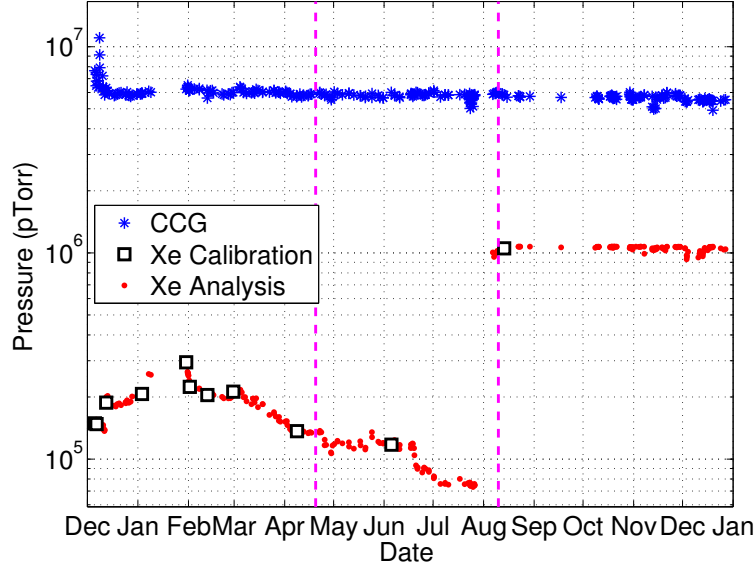


Figure 4: Xenon vapor pressure as measured by the RGA and CCG. The dashed, vertical magenta lines represent the course of the LUX underground science run used for the first WIMP limit published in [ref].

4.3. Correcting for Flow Rate

The measurements described in this article differ from that of Ref. [coldtrap paper] in that only a modest quantity of xenon gas is used for analysis in each sample, typically one standard liter, note that since xenon freezes in the coldtrap the sample is entirely recovered after the analysis. The sample pressure noticeably decreases during the measurement process, and also the mean flow rate varies between subsequent samples with the initial sample pressure. We account for the varying flow rate as follows. We normalize the average integral of the RGA signal to the integral of the flow rate, the leak rate normalized RGA partial pressure has been found to be constant in previous calibration [Kr paper, EXO ref]. Therefore the integrated RGA partial pressure of the impurity, which are measured in the first three minutes of flow, also depend on the average flow rate as the number of impurities that pass by the RGA per unit time is proportional to the flow. As shown in Figure 5, we have measured this dependence for gain corrected N_2 , He, Ar, CH_4 , and Kr signals, at our nominal

flow settings, using xenon gas from our calibration cylinder over a several month period. The data is well described by a simple linear fit for the dependence of gain corrected signal vs. mean flow rate.

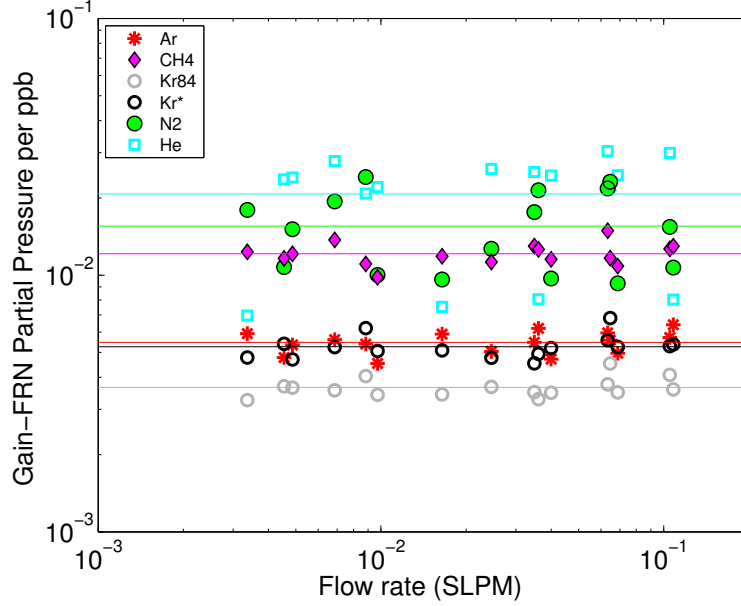


Figure 5: Gain corrected and leak rate normalized signals vs. flow rate for calibration xenon over eight months. The dominant source of variation is the RGA gain drift over time. (N2 and He calibration off by 2x, CH4, Ar, Kr are within 30%).

4.4. Gain and Flow Corrected Figure of Merit

To monitor for daily variations in the system response, including changes in the RGA gain, we periodically calibrate the analysis system using our cylinder of calibration xenon gas. After performing a periodic calibration we normalize the integral of the RGA partial pressure for a given species to both flow and the xenon partial pressure. First, the flow normalization is to remove the linear dependance of the signal on flow. Secondly, we correct for a drift in the RGA gain by normalizing to the partial pressure of xenon which is the one true constant between all samples as it represents the xenon ice vapor pressure at 77 K in the coldtrap system (given a fixed CCG pressure).

The purity measurement is made as follows, a full RGA scan of a typical calibration is shown in Figure 3. First, before the flow is induced the backgrounds of all impurity species are averaged for 60 seconds to record the baselines. Second, a command is sent to actuate a pneumatic valve opening a path from the sampling volume through one of the three leak valves into the coldtrap and RGA (see Figure []). Once the valve feedback confirms the open state, the code waits 30 seconds to

establish a stable flow and mitigate any transient effects. Third, all RGA signals and the flow rate are averaged over 150 seconds while the xenon flows from the sampling volume through the leak valve. During the measurement time about 60-70% of the sample is consumed and as the sample pressure depletes the flow rate diminishes thus, it is not practical to drag out the sample measurement to the last drop. Finally, the pneumatic value in front of the leak valve is shut cutting off flow to the cold trap. The signals on the RGA are measured for an additional ten minutes to ensure that the initial background levels return. It is important to note that all measurement timing is controlled with pneumatic valves via a C script on a slow control computer, also the flow settings dialed in on the three leak valves are permanently fixed. These measures were taken for the LUX analysis system to help mitigate systematic uncertainties associated with a human opening and closing the sensitive leak valves at measurement. In our early trials the flow of a given dial setting, after closing and reopening, was found to vary as much as 20% [coldtrap ref].

For a known initial mass fraction (ω_i [g/g]) the calibration figure of merit is defined by equation 2, defined as the gain corrected leak rate normalized RGA signal per $1 \cdot 10^{-9}$ g/g for a given impurity species. All subsequent purity measurements are normalized to flow rate and then corrected gain by taking the ratio of xenon vapor pressure. The purity measurement is divided by the previous calibration figure of merit and converted into a impurity concentration value in units of ppb ($1 \cdot 10^{-9}$ g/g), given in equation 3. We take the variation in calibration figure of merit before and after a sample for each individual impurity species as a systematic uncertainty in the measurement. This uncertainty represents systematics in gain drift, the RGA's response to xenon vapor pressure, the potential for small errors in calibration.

$$\omega_i = \frac{m_i[\text{g}]}{m_{\text{Xe}}[\text{g}]} \quad (1)$$

$$\text{Cal}_i = \langle \text{RGA}[\text{pTorr}] \rangle \cdot \langle \text{Flow}[\text{SLPM}] \rangle^{-1} \cdot \langle \text{XeP}[\text{pTorr}] \rangle^{-1} \cdot \langle \omega_i[\text{ppb}] \rangle^{-1} \quad (2)$$

$$\omega_j = \langle \text{RGA}[\text{pTorr}] \rangle \cdot \langle \text{Flow}[\text{SLPM}] \rangle^{-1} \cdot \langle \text{XeP}[\text{pTorr}] \rangle^{-1} \cdot \langle \text{Cal}_i[\text{ppb/SLPM}] \rangle^{-1} \quad (3)$$

Using Gain corrected and flow rate normalized calibration figure of merit (G-FRN) the measurements for a given impurity concentration are repeatable over a range of flows, RGA gains and over several months of calibrations. Figure 6 shows the variations in the system response to various impurities for our calibration xenon over a period of eight months with and without the xenon vapor pressure gain correction. The calibration figure of merit per ppb (g/g) for several impurities over a eight month period, with and without the additional gain normalization is also listed in Table 3. We find that the xenon vapor pressure measured by the CCG is stable to within 3% for all calibration,

as it should since it indirectly measures the vapor pressure of xenon ice at 77K. The xenon pressure, as measured by the RGA, noticeably drifts and tracks the other mass measurements while the CCG pressure remains fixed, indicating changes in the RGA gain. When normalizing the leak rate corrected figure of merit to the xenon vapor pressure there is a large improvement in resolution. The resolution relative to the mean improves from being over 100% to about 10%, with the exception of N₂ and He which may have additional systematics discussed in [link]. We find that over eight months of calibrations the calibration figure of merit per ppb for N₂, Ar, Kr, CH₄, He are stable to within 39%, 10%, 11%, 9%, 41% respectively.

Species	AMU	Cal _{FM} w/o gain μ	Cal _{FM} w/o gain σ/μ	Cal _{FM} $\mu \times 1000$	Cal _{FM} σ/μ
N ₂	28	5720	1.47	16.5	0.39
Ar	40	1840	1.08	6.0	0.10
Kr	84 (57.0%)	1290	1.20	3.9	0.11
Kr*	82,84,86 (85.6%)	1890	1.24	5.7	0.12
CH ₄	15	3890	1.01	13.0	0.09
He	4	6670	1.06	21.6	0.41
		$\langle \text{Pressure} \rangle$ (pTorr)	σ/μ		
Xe	132	3.04×10^5	1.01	–	–
CCG	n/a	6.03×10^6	0.031	–	–

Table 3: The calibration figure of merit per ppb (g/g) for several impurities over a eight month period, with and without the additional gain normalization (also shown in Figure 6). We find that the xenon vapor pressure measured by the CCG is stable to within 3% for all calibration. When normalizing the figure of merit to the xenon vapor pressure there is a large improvement in resolution. The resolution relative to the mean improves from being over 100% to about 10%, with the exception of N₂ and He which may have additional systematics discussed in [link].

4.5. O₂ and H₂ Figure of merit

Preparing xenon samples with low concentrations of O₂ and H₂ have proven to be difficult. Mainly because over time as the calibration bottle sits idle the stainless steel container absorbs O₂ and continuously out-gasses H₂. To account for the species-dependent response of the analysis system, including factors such as the probability that each impurity survives the cold trap and its ionization potential at the RGA filament, we apply a relative species response factor relative to nitrogen. For a given concentration of N₂ we have empirically observed that the equivalent RGA response is suppressed by a factor of 1.6 [ref] further, we have shows the linear response of the system to O₂ from $0.1 \cdot 10^{-9}$ g/g

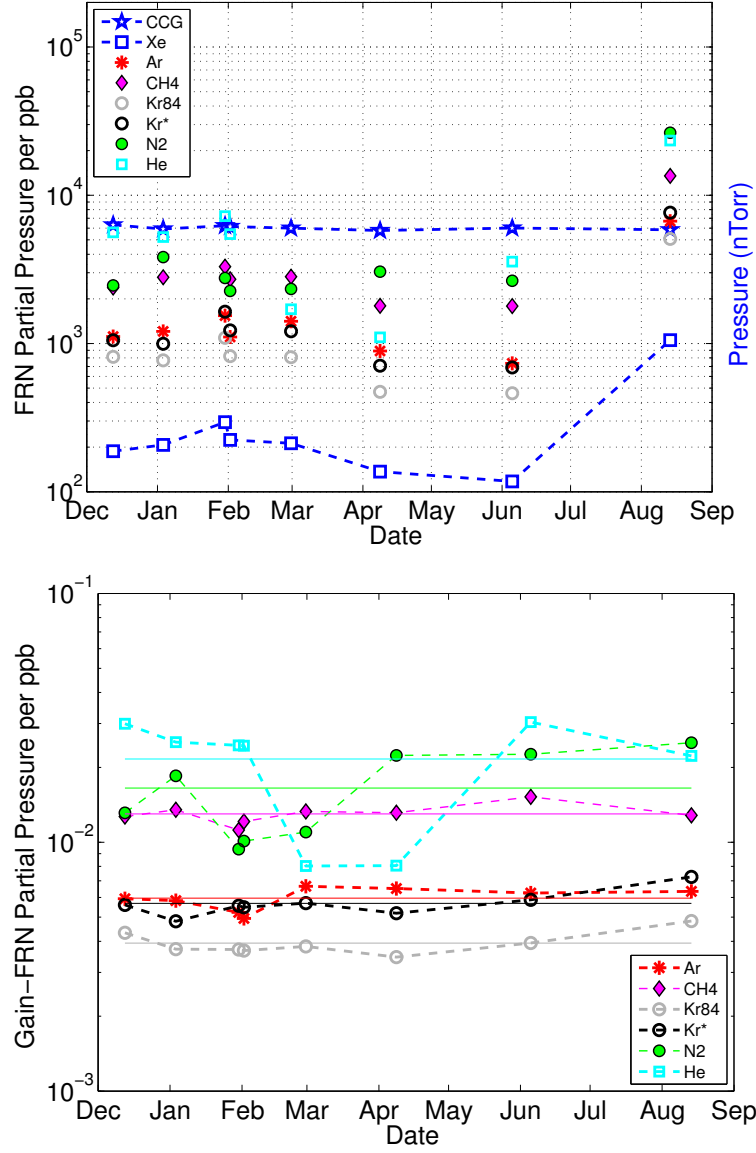


Figure 6: Top: Leak rate normalized RGA signals for several calibration species over eight months, the left axis shows the xenon vapor pressure measured by the RGA and cold cathode gauge (CCG). There is a strong correlation between the xenon vapor pressure and the leak rate normalized figure of merit. Bottom: Gain corrected and leak rate normalized figure of merit for calibration xenon over eight months.

to $1 \cdot 10^{-6}$ g/g [ref]. We have tested this conversion factor by calibrating with a sample prepared with $1 \cdot 10^{-6}$ g/g of N_2 and O_2 stored in an aluminum cylinder to minimize O_2 absorption. The application of this conversion factor from N_2 to O_2 response is valid to within 30% and has always been sensible for samples containing a slight known air contamination from which the N_2 , O_2 , Ar ratio can be used to scrutinize the calibration, constituting 75.5%, 23.2% and 1.3% of air by mass, respectively [ref]. (show measurements from Nov-Dec 2012 calibrations with O_2 and N_2 ?)

We also approximate the H_2 concentration of every sample multiplying the helium figure of merit by a by a factor, corresponding to the factor of two difference in the molar mass relative to helium when measuring in units of ppb g/g. Naively the RGA partial pressure response to a given mass concentration for a fixed ionization potential should be proportional to the number of moles of a given species. The hydrogen measurement may be off by up to a factor of four however, H_2 out-gasses indefinitely from stainless steel and it is good practice to track it's relative abundance overtime.

5. Applications of Gas Purity Analysis

We have found many useful applications for having the automated xenon sampling system integrated with the LUX circulation system. Now that we have discussed the methodology and calibrations we can discuss the results of the over 300 samples taken over the course of the underground science run. The primary purpose of the gas sampling system was to screen xenon for the accumulation of radioactive nobel isotopes present in the air, such as ^{85}Kr and ^{39}Ar for which dark matter experiments go through great lengths to remove from commercially available xenon [ref]. Prior to filling the LUX detector all xenon storage cylinders were screened to ensure they met purity and radioactive background requirements, then the xenon was monitored daily over the course of the science run. For the case of an air leak into the xenon circulation system impurities that degrade electron drift length, such as O_2 N_2 , are continuously removed by standard gettering technology [] while noble elements such as Ar and Kr accumulate. Experiments that only monitor liquid purity via the electron drift length are thus blind to small air leaks until the radioactive decay of ^{85}Kr begins to show up in data analysis, this scenario cost the Xenon100 experiment a year for Kr removal. The secondary purpose of the gas sampling system is to track the removal of electronegative impurities such as N_2 and O_2 that impede electron drift in the liquid xenon, a large drift length (~ 0.5 m) is necessary for TPC operation. The roll of the getter, installed in all xenon TPCs, is to remove electronegative impurities and hydrocarbons that outgas from detector components. The getters used can become saturated with impurities or suffer a loss of performance at high flow rate []. We have place sampling ports at the inlet and outlet of the getter and use weekly samples to monitor the getter's health by monitoring the penetration of N_2 , O_2 or hydrocarbons via H_2 and CH_4 . During a methane injection, prior to our tritiated methane calibrations, we monitored for the presence of H and C at the getters outlet to ensure the tritiated methane wouldn't be cracked by the getter and released as bare tritium into the LUX detector, bare tritium being 10x more penetrating than tritiated methane into plastics

[ref]. Finally, we have been able to use the sampling system to track the emanation of several species from detector components and measure solubility of impurities in xenon. Prior to the science run the detector was filled with gaseous xenon and circulated for a month, during which time extensive studies were done on the emanation of impurities for the teflon and polyethylene components. Solubilities for the other noble gases in liquid xenon could also be extracted by comparing the initial concentration of the xenon storage bottles to the concentration in the xenon gas after liquefaction. Limited data exists on solubilities of impurities in xenon and the measurements made with the LUX detector will be crucial for inferring liquid purity from the measured gas phase purity in future experiments.

Rn diffusion from air? ... we can set a limit from the Ar accumulation.

5.1. Xenon Inventory Screening

- ...plots for the Case gas chromatography
- ... Post Run03 compression and pre-Run04 results?
- measurements of bottles onsite.

5.2. Solubility of Impurities in Liquid Xenon

5.3. ^{85}Kr Monitoring

- ... plots for LUX science run

5.4. Gas Run Results, Outgassing